2 CONNECTIVITY CONFIGURATIONS


OPTION 1

DIRECT ACCESS between the units using an RJ45 cable. Preferred solution when **only one Eye Cubed system must be connected.**

OPTION 2

NETWORK ACCESS (Medical Network Access). Preferred solution when **several Eye Cubed units must be connected.**

QUICK START: EYE CUBED VIEWER ON ABSolu UNIT

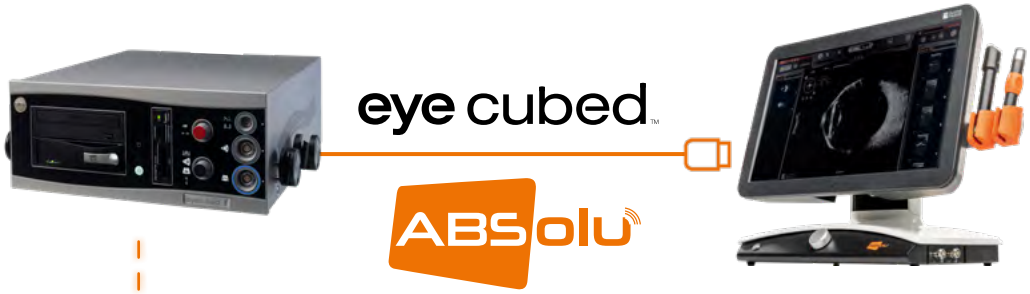

FOLLOW THE GUIDELINE

1. Switch on the Eye Cubed and ABSolu units


2. Click on the key of ABSolu keyboard to access to the task bar

3. Click on the Remote Desktop Connection icon located in the task bar to launch the “Remote Access Connection” software

   The “Remote Access Connection” function is available on Windows and is not a function specific to the ABSolu software. It’s easy to implement (to be activated on ABSolu and on Eye Cubed) – “Remote Access Connection” icon located in the taskbar to switch from one ultrasound platform to another one.

4. Select the needed computer in the list:

5. Click on the “Connect” key, the following screen will be displayed:

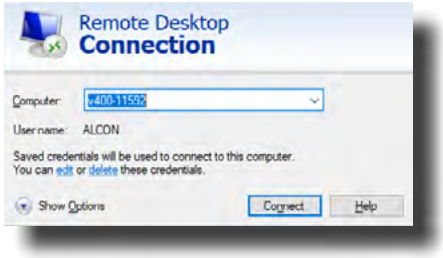

6. Fill the password field (if a password is set in the Eye Cubed Computer, if not click directly on “OK” key):

   In order to avoid to fill in the password at each connection, click on the “Remember me”. This selection will allow you that the above window will no more displayed.

7. Click on OK and the following screen will be displayed:

   The system may display the following error message: “User name or password are incorrect” (due to empty password), click on the “OK” key.



   To avoid this error message, a password is needed in the Windows Session of the Eye Cubed computer: refer to the Frequently Asked Questions chapter.

8. The ABSolu unit is connected to the Eye Cubed computer and the data of the Eye Cubed computer screen are displayed on the ABSolu screen.

9. Use the Tab + Alt keys to display Eye Cubed computer screen. Minimize the “Remote Desktop Connection” window on the menu barre at the top of the screen to display ABSolu screen.

   Check that the Remote Desktop Access function is not displayed in full screen and that the “Smart sizing” selection is selected in the remote connection to display all functionalities of the Eye Cubed unit.

A slight delay may happen when viewing the Eye Cubed software on the ABSolu screen. Do not perform scans when using the ABSolu screen as a viewer.

Quantel medical BY LUMIBIRD MEDICAL

> For more detailed information about this procedure, please refer to the technical notice to install Eye Cubed Viewer on ABSolu unit.