2 CONNECTIVITY CONFIGURATIONS

OPTION 1

DIRECT ACCESS between the units using an RJ45 cable. Preferred solution when **only one Aviso system must be connected.**

OPTION 2

NETWORK ACCESS (Medical Network Access). Preferred solution when **several Aviso units could be connected.**

QUICK START: AVISO® VIEWER ON ABSolu® UNIT



FOLLOW THE GUIDELINE


1. Switch on the Aviso and ABSolu units

2. Click on the key of ABSolu keyboard to access to the task bar

3. Click on the Remote Desktop Connection icon located in the task bar to launch the “Remote Access Connection” software

The “Remote Access Connection” function is available on Windows and is not a function specific to the ABSolu software. It’s easy to implement (to be activated on ABSolu and on Aviso) – “Remote Access Connection” icon located in the taskbar to switch from one ultrasound platform to another one.

4. Select the needed computer in the list:

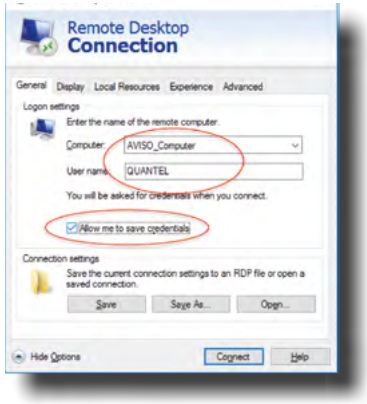


5. Click on the “Connect” key, the following screen will be displayed:

6. Fill the password field (if a password is set in the Aviso Computer, if not click directly on “OK” key):

In order to avoid to fill in the password at each connection, click on the “Remember me”. This selection will allow you that the above window will no more displayed.

7. Click on OK and the following screen will be displayed:

The system may display the following error message: “User name or password are incorrect” (due to empty password), click on the “OK” key.

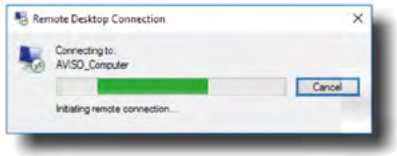

To avoid this error message, a password is needed in the Windows Session of the Aviso computer: refer to the Frequently Asked Questions chapter.

8. The ABSolu unit is connected to the Aviso computer and the data of the Aviso computer screen are displayed on the ABSolu screen.

9. Use the Tab + Alt keys to display Aviso computer screen. Minimize the “Remote Desktop Connection” window on the menu barre at the top of the screen to display ABSolu screen.

Check that the Remote Desktop Access function is not displayed in full screen and that the “Smart sizing” selection is selected in the remote connection to display all functionalities of the Aviso unit.

A slight delay may happen when viewing the Aviso software on the ABSolu screen. Do not perform scans when using the ABSolu screen as a viewer.

Quantel medical BY LUMIBIRD MEDICAL

> For more detailed information about this procedure, please refer to the technical notice to install Aviso Viewer on ABSolu unit.